


Air Ankle Stirrups


Application in four easy steps

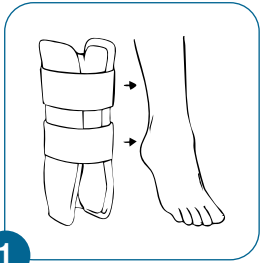

1

Position the Stirrup


Place the stirrup around the ankle, rigid shells on both sides.

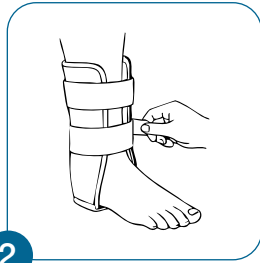

2

Secure the Straps


Fasten the hook-and-loop straps securely around the ankle for a comfortable fit.

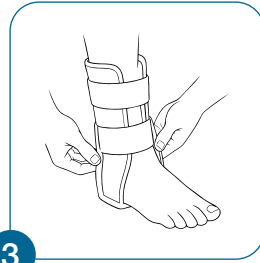

3

Adjust Air Support

Adjust the inflatable air bladders to provide a secure, customized fit.

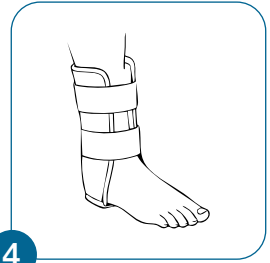

4

Check Support

Ensure the stirrup is comfortably positioned and provides ankle support.

Reusable Design:

Repeated use. Applies to Medzah Air Ankle Stirrup, Item Number: 4000-00.

Contents:

Polyester fabric, nylon webbing, polypropylene, polyurethane foam, TPU bladder, nylon, polyester, and plastic components.

Caution:

Inspect the stirrup, straps, and air bladders before use. Do not use if damaged. Ensure the device is securely fitted without excessive compression or discomfort. Follow a healthcare professional's directions.

Storage and Handling:

Store at room temperature in a clean, dry location. Keep away from excessive heat, moisture, and sharp objects. Inspect the device before each use.

Important Use Information:

Inspect the stirrup, straps, and air bladders before use. Position the device securely around the ankle and adjust the hook-and-loop straps and air bladders for a comfortable fit. Ensure proper support without restricting circulation.

