


Alcohol Prep Pads


Application in three easy steps

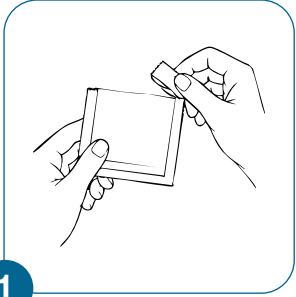

1

Open the Pad

Tear open the individual wrapper and remove the sterile alcohol prep pad.

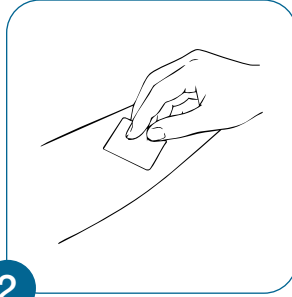

2

Prepare the Skin

Wipe the area thoroughly with the pad and allow the skin to dry before the procedure.



3

Discard After Use

Dispose of the used pad immediately after use for hygiene.

Single-Use Design:

Single use. Applies to Medzah Alcohol Prep Pads, Item Number: ALC001 (Medium), ALC002 (Large).

Contents:

Non-woven fabric saturated with 70% v/v isopropyl alcohol.

Caution:

For external use only. Do not use if the package is opened or damaged. Avoid contact with eyes. Do not use on broken skin. Do not reuse the pad. Keep out of reach of children. Flammable. Keep away from fire or flame.

Storage and Handling:

Store in a clean, dry location at room temperature. Keep individually wrapped pads sealed until use and protect from excessive heat and moisture.

Important Use Information:

Inspect the individual wrapper before use. Do not use if the package is opened or damaged. Use the sterile pad to prepare the skin before injections, blood draws, IV insertions, and other medical procedures.

