


Aluminum Quad Canes


Use in three easy steps

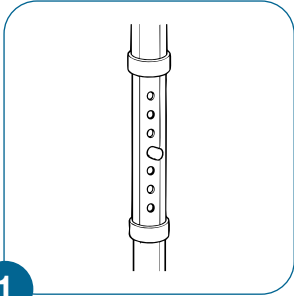

1

Adjust the Height


Set the cane to a comfortable height between 27–37 inches for proper posture.

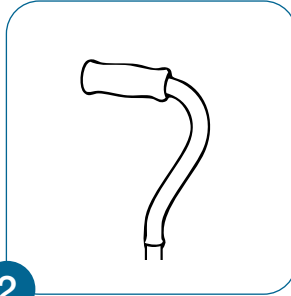

2

Check the Grip

Ensure the PVC handle provides a secure, comfortable grip.

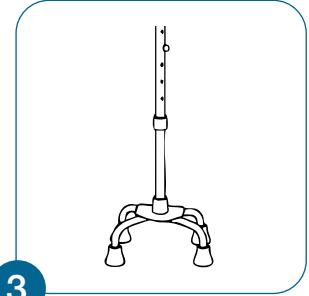

3

Walk with Support

Hold the handle firmly and place the quad base securely on the ground while walking.

Reusable Design:

Repeated uses. Applies to Medzah Aluminum Quad Canes, Item Number: MUJI1210-01 (Small Base), MUJI1210-02 (Large Base).

Contents:

Aluminum quad cane with PVC offset handle and quad base.

Caution:

Do not use if the cane is damaged or if the height adjustment is not securely locked. Ensure the cane is properly adjusted to provide stable support. Follow a healthcare professional's directions. Maximum weight capacity: 300 lbs. Do not exceed.

Storage and Handling:

Store in a clean, dry location at room temperature. Handle carefully to prevent damage to the cane or handle. Inspect the cane before each use.

Important Use Information:

Inspect the cane before use. Adjust the height to provide proper support and ensure the adjustment is securely locked. Use the handle for comfortable handling and follow a healthcare professional's directions.

