


Bionex Nitrile Gloves


Application in three easy steps

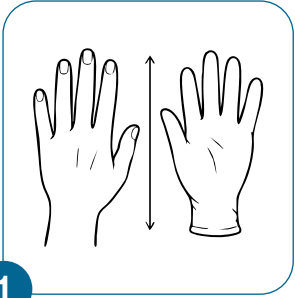

1

Select the Right Size

Choose a glove size that fits your hand snugly and comfortably.

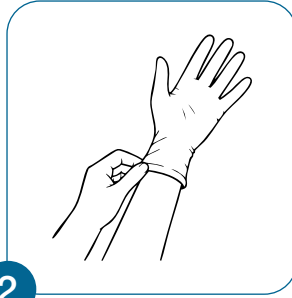

2

Don the Glove

Gently pull the glove on and fit it smoothly over the fingers, palm, and wrist.



3

Remove and Dispose

Carefully peel the glove off inside out and discard it after a single use.

Single-Use Design:

Single use. Applies to Medzah Bionex Nitrile Gloves, Item Numbers: BNMD01 (XS), BNMD02 (S), BNMD03 (M), BNMD04 (L), BNMD05 (XL).

Contents:

Nitrile examination glove. Length 9.5 in; finger thickness 0.10 mm (4.0 mil); palm thickness 0.07 mm (3.0 mil).

Caution:

Single use only. Do not reuse. Do not use if the glove or package is damaged. Discontinue use if skin irritation occurs.

Storage and Handling:

Store Bionex Nitrile Gloves in a clean, dry location at room temperature. Keep away from excessive moisture, heat, direct sunlight, and sharp objects that may damage the product or packaging.

Important Use Information:

Select the appropriate size for a snug, comfortable fit. Dispense one glove at a time, cuff-first. Tested for use with chemotherapy drugs per ASTM D6978-05. Some chemotherapy drugs have short breakthrough times; refer to the spec sheet before use. Dispose of used gloves according to local regulations.

