


Cervical Collar


Application in three easy steps

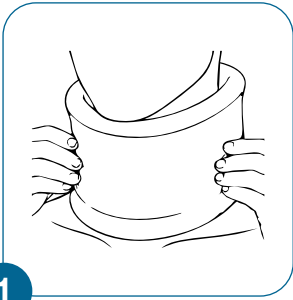

1

Position the Collar


Place the collar comfortably around the neck with the chin supported.

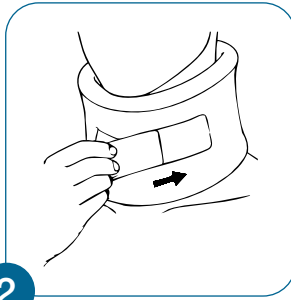

2

Secure the Closure

Fasten the hook-and-loop closure securely for a comfortable fit.

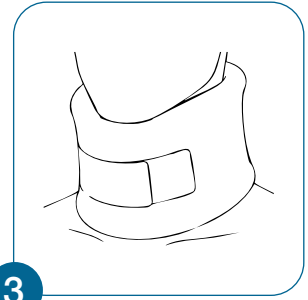

3

Check the Support

Ensure the collar provides comfortable support while reducing neck movement.

Reusable Design:

Repeated use. Applies to Medzah Cervical Collar, Item Number: 6000-25, 6000-30, 6000-40.

Contents:

Open-cell polyurethane foam, polyester and spandex jersey knit fabric, polyester hook-and-loop, and polyester straps.

Caution:

Do not use if the collar or closure is damaged. Ensure the collar is securely fitted without excessive pressure or discomfort. Follow a healthcare professional's directions.

Storage and Handling:

Store at room temperature in a clean, dry location. Keep away from excessive moisture, heat, and sharp objects. Inspect the collar before each use.

Important Use Information:

Inspect the collar and closure before use. Position the collar comfortably around the neck and secure the hook-and-loop closure for a custom fit. Ensure the collar provides support without causing excessive pressure or discomfort.

