


Compression Knee Immobilizer


Application in three easy steps

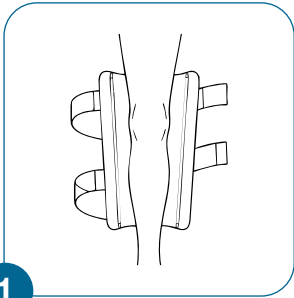

1

Position the Immobilizer


Secure Immobilizer around the knee, reinforced stays along both sides of the leg.

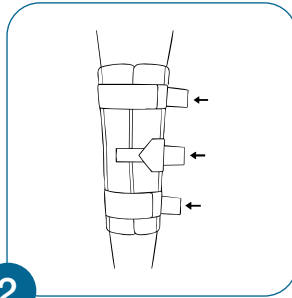

2

Secure the Straps

Fasten the adjustable elastic straps to provide compression and maintain secure positioning.

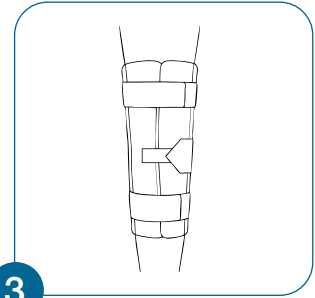

3

Check the Fit

Ensure the immobilizer is snug, the stays are aligned, and circulation is not restricted.

Reusable Design:

Repeated use. Applies to Medzah Compression Knee Immobilizer, Item Numbers: 2000-00 (12 in), 2000-00 (16 in), 2000-00 (20 in), 2000-00 (24 in).

Contents:

Open-cell polyurethane foam, polyester-spandex jersey knit fabric, and aluminum stays.

Caution:

Do not use if the immobilizer, stays, or straps are damaged. Ensure the immobilizer is securely fitted without excessive pressure or discomfort. Follow a healthcare professional's directions.

Storage and Handling:

Store at room temperature in a clean, dry location. Keep away from excessive moisture, heat, and sharp objects. Inspect the immobilizer before each use.

Important Use Information:

Position the contoured reinforced stays along the knee and adjust the elastic straps to provide even compression and secure support. The material may be trimmed to accommodate individual requirements.

