


Elbow Splint


Application in three easy steps

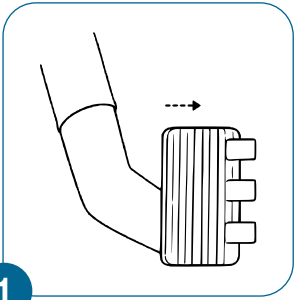

1

Position the Splint


Place the splint around the elbow and align it to protect the IV site.

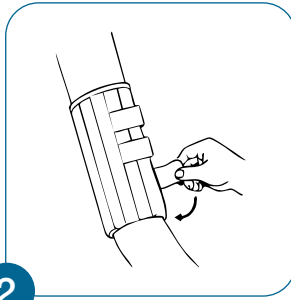

2

Adjust the Splint

Fasten the hook-and-loop tabs securely around the arm.

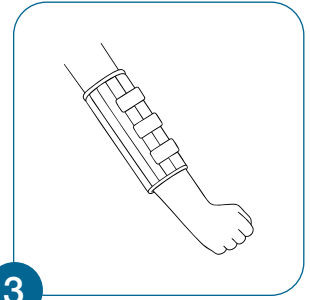

3

Check Placement

Ensure the IV site remains protected while allowing arm movement.

Reusable Design:

Repeated use. Applies to Medzah Elbow Splint, Item Number: MED-1001, MED-1002, MED-1003, MED-1004.

Contents:

Nylon-laminated foam with reinforced aluminum supports and hook-and-loop tab closures.

Caution:

Do not use if the splint is damaged. Ensure the splint is securely positioned without restricting circulation. Check the IV site regularly and adjust or remove the splint if discomfort occurs. Follow a healthcare professional's directions.

Storage and Handling:

Store in a clean, dry location at room temperature. Keep away from excessive moisture, heat, and sharp objects. Inspect the splint before each use.

Important Use Information:

Inspect the splint and closures before use. Position securely around the elbow to protect the IV site while allowing arm movement. Ensure the fit remains comfortable and does not interfere with circulation.

