


Infant Heel Warmer


Activation in four easy steps

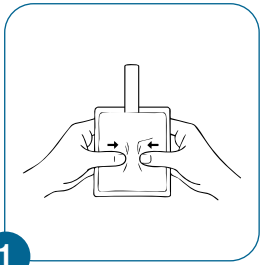

1

Squeeze the Pack


Squeeze the pack to activate the inner fluid pouch.

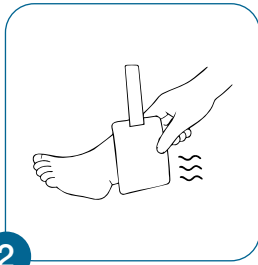

2

Apply Warmer Therapy


Place the warmer over the infant's heel for localized warmth.

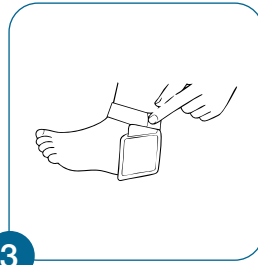

3

Secure the Warmer

Use the attached tape to secure the warmer firmly in place.

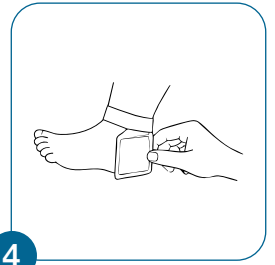

4

Check Fit & Monitor

Ensure the warmer is comfortable and monitor the infant.

Disposal after single use:

Single use. Applies to Medzah Infant Heel Warmer, Item Number: SAM-108.

Contents:

Food-grade sodium acetate and water.

Caution:

Do not puncture or damage the warmer. Do not use if damaged or leaking. For external use only. Monitor the infant's skin during use for discomfort or adverse reactions. Discontinue use if any reaction occurs.

Storage and Handling:

Store at room temperature. Keep away from extreme temperatures.

Important Use Information:

Inspect the warmer before use. Do not use if damaged or leaking. Activate by squeezing the inner fluid pouch and allow the warmer to reach 101°F–104°F (38°C–40°C) before application. Follow a healthcare professional's directions.

