



# Instant Cold Pack

## Activation in four easy steps

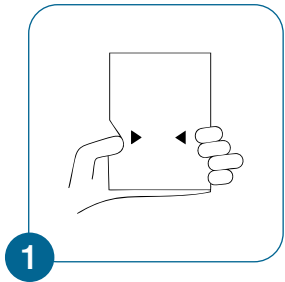


### 1 Squeeze the Pack

Squeeze the cold pack firmly to activate the internal barrier.

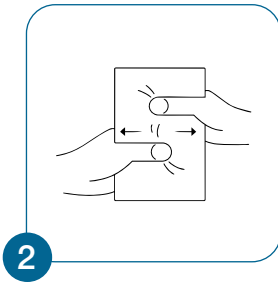


### 2 Shake the Pack

Shake the activated pack thoroughly to distribute the contents evenly.

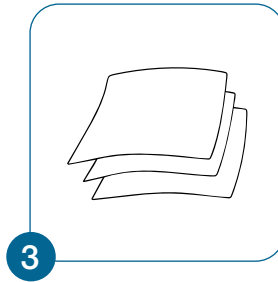

### 3 Wrap in Cloth

Wrap the activated pack in a cloth before applying the pack.

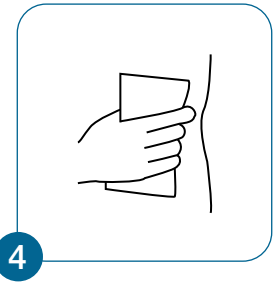

### 4 Apply Cold Therapy

Apply the wrapped pack to the affected area for localized cold therapy.

### Disposal after single use:

Single use. Applies to Medzah Instant Cold Pack, Item Number: SAM-1001.

### Contents:

Calcium ammonium nitrate and water.

### Caution:

Do not use on infants. Do not puncture or strike the pack in any manner. Do not use if punctured, torn, or leaking. For external use only. Do not apply directly to the skin. Monitor the skin for discomfort or adverse reactions during use. Discontinue use if any reaction occurs.

### Storage and Handling:

Store at room temperature. Do not refrigerate or freeze prior to use. Keep away from extreme temperature conditions. Follow the activation instructions printed on the pack before application.

### Important Use Information:

Inspect the pack before use. Do not use if torn, punctured, leaking, or damaged. Wrap the activated pack in a cloth before application and follow a healthcare professional's directions. Provides approximately 20 minutes of cold therapy.

