


Instant Hot Pack


Activation in three easy steps

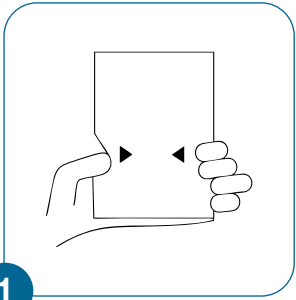

1

Squeeze the Pack


Squeeze the hot pack firmly to activate the internal barrier.

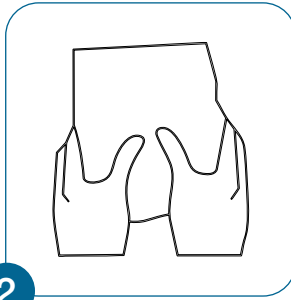

2

Knead the Pack

Knead the pack thoroughly to activate the warming reaction.

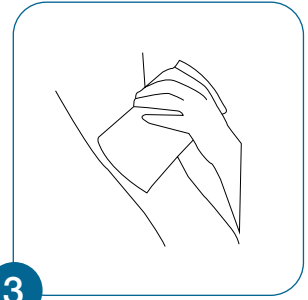

3

Apply Heat Therapy

Apply the pack with a protective layer or towel between the pack and skin.

Disposal after single use:

Single use. Applies to Medzah Instant Hot Pack, Item Number: SAM-2000.

Contents:

Sodium acetate and water.

Caution:

Do not puncture or damage the pack. Do not use if torn, punctured, leaking, or damaged. For external use only. Use with a protective layer or towel between the pack and skin. Monitor the skin during use for discomfort or adverse reactions. Discontinue use if any reaction occurs.

Storage and Handling:

Store at room temperature. Keep away from extreme temperatures. Do not use after the expiration date (3-year shelf life).

Important Use Information:

Inspect the pack before use. Do not use if torn, punctured, leaking, or damaged. Knead the pack to activate. Apply with a protective layer or towel between the pack and skin, and follow a healthcare professional's directions. Provides approximately 20 minutes of heat at 105°F–110°F.

