


Medmax Underpads


Application in three easy steps

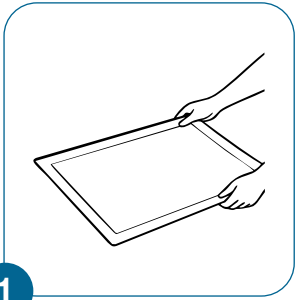

1

Position the Underpad


Place the underpad flat with the blue indicator line facing down.

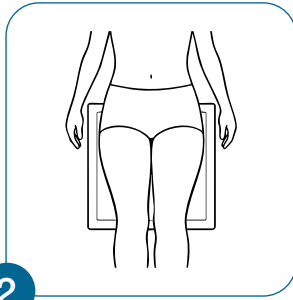

2

Align the Pad

Center the pad under the patient over the area to be protected.

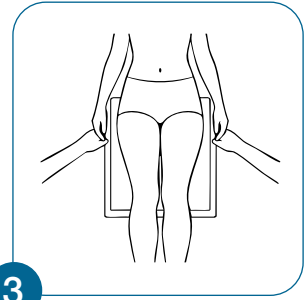

3

Smooth and Check

Smooth out any folds so the pad lies flat and stays in place during movement.

Single-Use Design:

Single use. Applies to Medzah Medmax Underpads, Item Numbers: IPUP02 (23 x 36 in), IPUP03 (30 x 36 in), IPUP04 (26 x 36 in), IPUP05 (30 x 36 in).

Contents:

Multi-layer underpad with hydrophilic non-woven fabric top layer, ultra soft tissue filling, fluff and super absorbent polymer core, and waterproof breathable PE film bottom layer. Latex-free.

Caution:

For external use only. Do not reuse. Do not use if the package or pad is damaged. Discontinue use if skin irritation occurs.

Storage and Handling:

Store Medmax Underpads in a clean, dry location at room temperature. Keep away from excessive moisture, heat, direct sunlight, and sharp objects that may damage the product or packaging.

Important Use Information:

Select the appropriate size according to the area to be protected. Place the pad with the blue indicator line facing down, smooth it flat, and ensure the pad remains in position during use. Weight capacity up to 400 lbs (IPUP03, IPUP05).

