


Non-Woven Instant Cold Pack


Activation in three easy steps

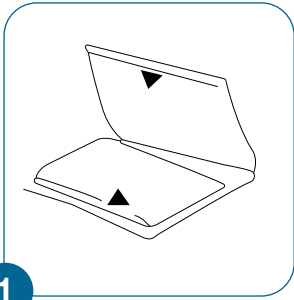

1

Fold the Pack


Fold the cold pack firmly in half at the indicated activation area.

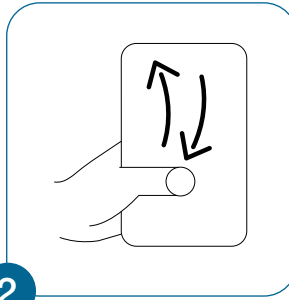

2

Shake the Pack

Shake the activated pack thoroughly to help distribute the cooling contents evenly.

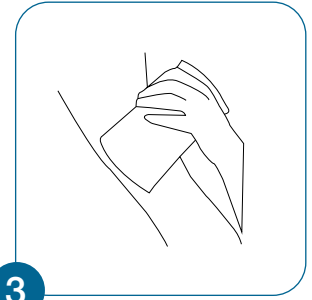

3

Apply Cold Therapy

Place the activated cold pack over the desired area for localized cold therapy.

Disposal after single use:

Single use. Applies to Medzah Non-Woven Instant Cold Pack, Item Number: SAM-3001.

Contents:

Calcium ammonium nitrate and water.

Caution:

Do not puncture or strike the pack in any manner. Do not use if punctured, torn, or leaking. For external use only. Individuals with circulatory problems should use only at the direction of a physician. Do not apply to broken skin. Monitor sensitive skin for discomfort, burns, or adverse reactions during use. Discontinue use if any of these reactions occur.

Storage and Handling:

Store at room temperature. Do not refrigerate or freeze prior to use. Do not heat. Keep away from extreme temperature conditions. Follow the activation instructions printed on the pack before application.

Important Use Information:

Inspect the pack before use. Do not use if torn, punctured, leaking, or damaged. Position comfortably and follow a healthcare professional's directions. Provides approximately 20 minutes of cold therapy and does not drop below 50°F.

