


Non-Woven Instant Hot Pack


Activation in three easy steps

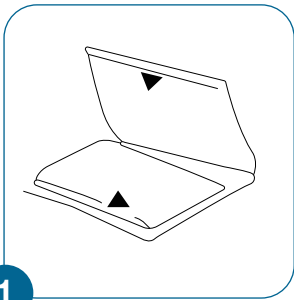

1

Fold the Pack


Fold the hot pack firmly in half at the indicated activation area.

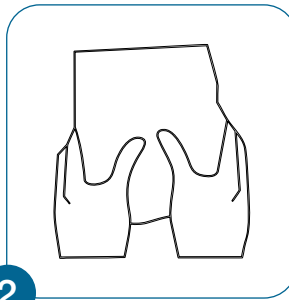

2

Knead the Pack

Knead the folded pack thoroughly to activate the warming reaction.

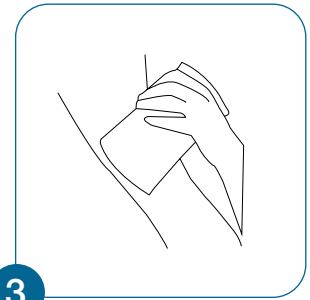

3

Apply Heat Therapy

Place the activated hot pack over the desired area for localized heat therapy.

Disposal after single use:

Single use. Applies to Medzah Non-Woven Instant Hot Pack, Item Number: SAM-2001.

Contents:

Sodium acetate and water.

Caution:

Do not puncture or damage the pack. Do not use if torn, punctured, leaking, or damaged. For external use only. Monitor the skin during use for discomfort or adverse reactions. Discontinue use if any reaction occurs.

Storage and Handling:

Store at room temperature. Keep away from extreme temperatures. Do not use after the expiration date (3-year shelf life).

Important Use Information:

Inspect the pack before use. Do not use if torn, punctured, leaking, or damaged. Fold and knead the pack to activate, then apply directly to the skin for localized heat therapy. Follow a healthcare professional's directions. Provides approximately 20 minutes of heat at 105°F–110°F.

