



Offset Handle Fashion Canes

Use in three easy steps

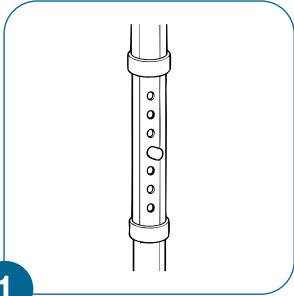


Adjust the Height

Set the cane to a comfortable height between 27–37 inches for proper posture.

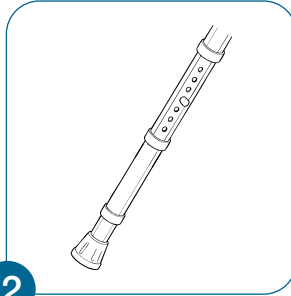

Secure the Adjustment

Check that the height adjustment is fully secured before use.

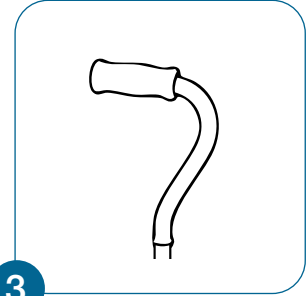

Grip the Handle

Hold the handle firmly and place the cane tip securely on the ground while walking.

Reusable Design:

Repeated uses. Applies to Medzah Offset Handle Fashion Cane, Item No. MUJI3210-0106.

Contents:

Aluminum cane with polypropylene handle and wrist strap.

Caution:

Do not use if the cane is damaged or if the height adjustment is not securely locked. Ensure the cane is properly adjusted and provides stable support. Follow a healthcare professional's directions. Maximum weight capacity: 220 lbs. Do not exceed.

Storage and Handling:

Store in a clean, dry location at room temperature. Handle carefully to prevent damage to the cane or handle. Inspect the cane before each use.

Important Use Information:

Inspect the cane before use. Adjust the height to provide a comfortable fit and ensure the adjustment is securely locked. Use the handle and wrist strap for comfortable handling and support.

