


Plastic Specimen Containers

Application in three easy steps



Collect the Specimen

Place the required specimen into the container without exceeding 90 mL capacity.

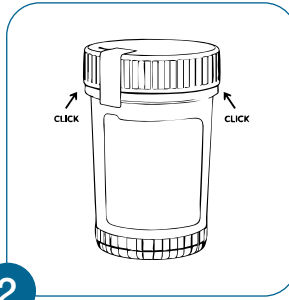

Secure the Cap

Press the cap firmly until both locking points click into place.



Check the Seal

Confirm the cap is locked and the tamper-evident seal is intact before transport.

Single-Use Design:

Single use. Applies to Medzah Plastic Specimen Containers, Item No. EVON101 (no temperature strip), EVON102 (with temperature strip).

Contents:

Medical-grade polypropylene container with double-click locking cap and integrated tamper-evident seal. EO sterilized. 95 kPa certified.

Caution:

Do not use if the container or cap is damaged. Do not exceed the 90 mL capacity. Ensure the cap is securely locked before handling or transport. Follow applicable specimen collection and handling procedures.

Storage and Handling:

Store in a clean, dry location at room temperature. Handle the container and cap carefully to prevent damage. Inspect the container for any signs of damage before use and ensure the cap is intact.

Important Use Information:

Inspect the container and cap before use. Collect the required specimen, secure the double-click locking cap, and confirm the seal is intact before transport. Label the container using the pre-printed panel.

