

Post-Op Shoe

Application in three easy steps



1

Open the Shoe


Open the hook-and-loop straps to allow easy foot placement.



2

Position the Foot

Place the foot comfortably inside the shoe with the heel positioned securely.

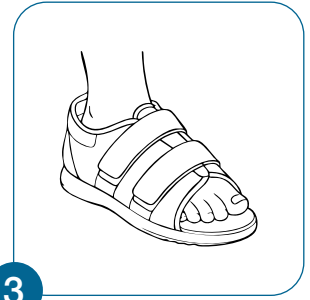

3

Secure the Shoe

Fasten the straps firmly for a comfortable and secure fit.

Reusable Design:

Repeated use. Applies to Medzah Post-Op Shoe, Item Number: 5000-01, 5000-02, 5000-03, 5000-04 (Male); 5010-01, 5010-02, 5010-03 (Female).

Contents:

Polyester mesh upper with PE foam and cotton, TPU heel counter, rigid rocker sole, and hook-and-loop straps.

Caution:

Do not use if the shoe or closure is damaged. Ensure the shoe is securely fitted without excessive pressure or discomfort. Follow a healthcare professional's directions.

Storage and Handling:

Store at room temperature in a clean, dry location. Keep away from excessive moisture, heat, and sharp objects. Inspect the shoe before each use.

Important Use Information:

Inspect the shoe and closure before use. Ensure the foot is positioned comfortably with the heel securely supported. Fasten the closure for a firm, comfortable fit and follow a healthcare professional's directions.

