


Shoulder Immobilizer


Application in four easy steps

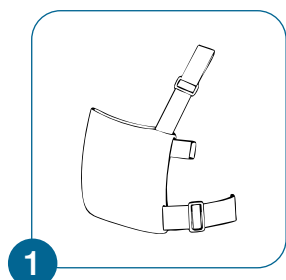

1

Prepare Immobilizer


Open the hook-and-loop closures and use the straps for easy placement.

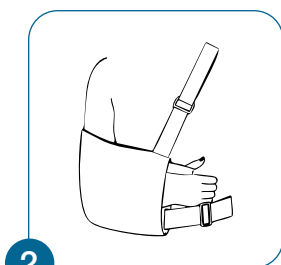

2

Place the Arm


Place the injured arm inside and adjust the hand position securely.

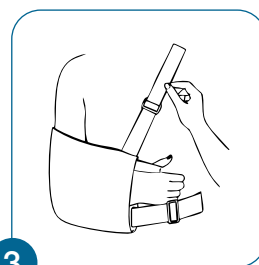

3

Adjust the Straps

Adjust straps for a personalized fit and position for comfort.

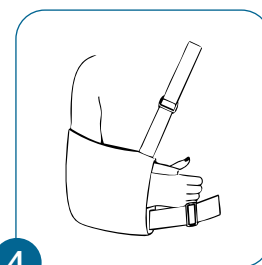

4

Secure Waist Strap

Fasten the waist strap to help prevent rotation and maintain immobilization.

Reusable Design:

Repeated use. Applies to Medzah Shoulder Immobilizer, Item Number: 1000-00, 1000-01, 1000-02, 1000-03, 1000-04.

Contents:

Cotton sling, polypropylene loops, polyester straps, foam straps, shoulder pad, and hook-and-loop closure components.

Caution:

Do not use if the immobilizer, sling, straps, or closures are damaged. Ensure the immobilizer is securely fitted without excessive pressure or discomfort. Follow a healthcare professional's directions.

Storage and Handling:

Store at room temperature in a clean, dry location. Keep away from excessive moisture, heat, and sharp objects. Inspect the immobilizer and straps before each use.

Important Use Information:

Position the arm comfortably in the sling and adjust all hook-and-loop closures for a secure, personalized fit. Use the waist strap to help limit shoulder rotation and maintain proper immobilization. The immobilizer is designed to fit either the right or left arm.

