


Splash Mask Face Shields


Use in three easy steps

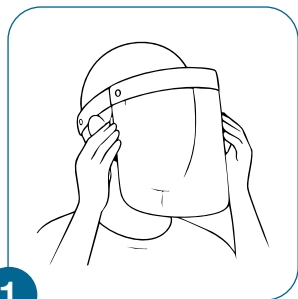

1

Position the Shield


Place the face shield in front of the face and ensure it provides full-face coverage.

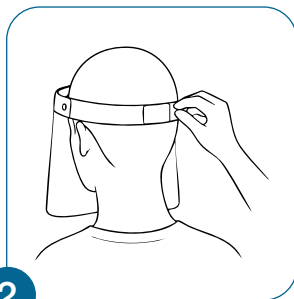

2

Secure the Strap

Fasten the elastic strap securely for a comfortable, stable fit.

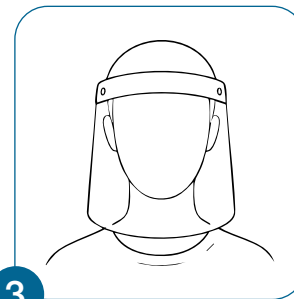

3

Check Visibility

Confirm the shield is properly positioned and provides clear visibility before use.

Reusable Design:

Repeated uses. Applies to Medzah Splash Mask Face Shields, Item Number: 11000-001 (Vented), 11000-003 (Solid).

Contents:

Splash Mask Face Shield made from latex-free, fiberglass-free materials.

Caution:

Do not use if the face shield or packaging is damaged or contaminated. Do not use if the face shield does not provide adequate coverage or visibility.

Storage and Handling:

Store in a clean, dry location at room temperature, away from direct sunlight, excessive heat, and moisture. Keep in its original packaging or a clean, protected area when not in use.

Important Use Information:

Inspect the face shield and packaging before use. Do not use if damaged or contaminated. Ensure the face shield is properly positioned and provides intended full-face coverage and clear visibility during use.

