


Steel Folding Walker With 5" Wheels


Use in three easy steps

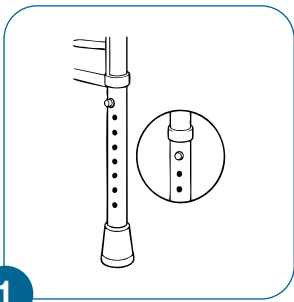

1

Adjust the Height


Set the walker to the desired height (32.5-38.5 in) using the push-button height settings.

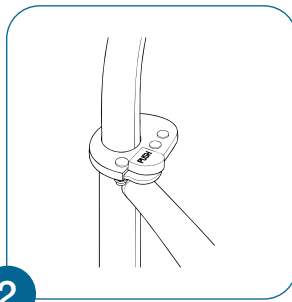

2

Open and Lock

Push the folding mechanism to open the walker and confirm it locks with an audible click.

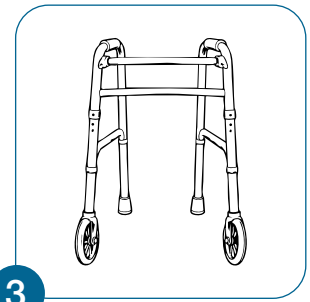

3

Check Stability

Ensure the walker is fully opened and stable before use.

Reusable Design:

Repeated uses. Applies to Medzah Steel Folding Walker With 5" Wheels, Item No. MDZH904SE, MDZH904S.

Contents:

Steel frame with PVC handles, 5" rubber wheels, and anti-skid rubber tips.

Caution:

Do not use if the walker is damaged or if the locking mechanism does not engage properly. Ensure all height adjustments are securely locked before use. Use on stable surfaces and follow a healthcare professional's directions. Maximum weight capacity: 400 lbs. Do not exceed.

Storage and Handling:

Store in a clean, dry location. Fold flat for transport and storage. Inspect the frame, wheels, handles, and locking mechanism before each use.

Important Use Information:

Inspect the walker before use. Adjust the height for a comfortable fit and ensure the locking mechanism is fully engaged. Confirm the walker is stable before walking and use the handles for support.

