


Thumb Spica


Application in three easy steps

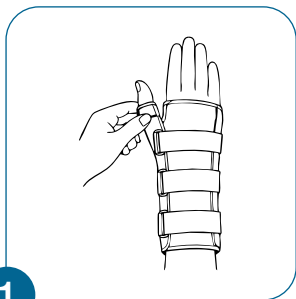

1

Position the Spica


Place the thumb in the neutral position with the hand comfortably supported.

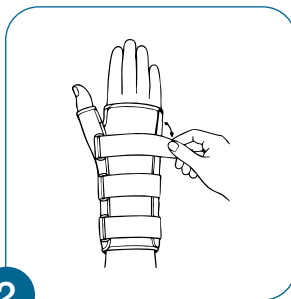

2

Adjust the Fit

Allow the expandable panel to conform to the patient's anatomy.

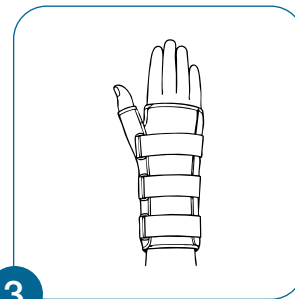

3

Secure the Strap

Fasten the adjustable auto-lock strap to secure the selected fit.

Reusable Design:

Repeated use. Applies to Medzah Thumb Spica, Item Number: 3010-10, 3010-11, 3010-12, 3010-20, 3010-21, 3010-22.

Contents:

Open-cell polyurethane foam with aluminum and nylon support stays, polyester, nylon, steel, and POM components.

Caution:

Do not use if the thumb spica, support stays, or strap are damaged. Ensure the thumb is positioned comfortably without excessive pressure or restriction of circulation. Follow a healthcare professional's directions.

Storage and Handling:

Store at room temperature in a clean, dry location. Keep away from excessive moisture, heat, and sharp objects. Inspect the thumb spica before each use.

Important Use Information:

Inspect the thumb spica and strap before use. Position the thumb in a neutral position and adjust the expandable panel for a custom fit. Secure the auto-lock strap and ensure the hand remains comfortably supported with functional movement.

Select the correct size and side (left or right) before use.

