


Tri-Panel Knee Immobilizer


Application in three easy steps

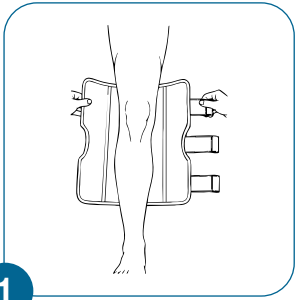

1

Position the Immobilizer


Place the immobilizer around the leg with the knee centered in the opening.

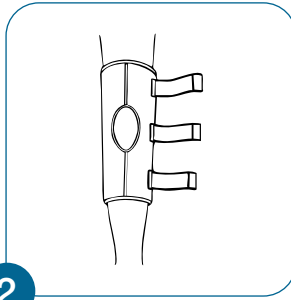

2

Adjust the Panels

Adjust the three panels to fit comfortably around the leg and accommodate swelling.

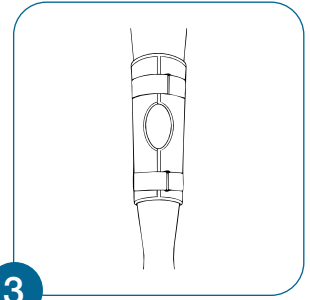

3

Secure the Straps

Fasten the elastic straps securely for a snug, supportive fit.

Reusable Design:

Repeated use. Applies to Medzah Tri-Panel Knee Immobilizer, Item Number: 40012, 40016, 40020, 40024.

Contents:

Tri-panel knee immobilizer with rigid/semi-rigid support components, soft padding, and elastic straps.

Caution:

Do not use if the immobilizer or straps are damaged. Ensure the knee is properly positioned and the straps are secured without excessive pressure or discomfort. Follow a healthcare professional's directions.

Storage and Handling:

Store at room temperature in a clean, dry location. Wash to maintain cleanliness and performance. Inspect the immobilizer before each use.

Important Use Information:

Inspect the immobilizer and straps before use. Position the knee centrally and adjust the panels to provide comfortable support. Secure the elastic straps evenly and ensure proper circulation during use.

