

Urine Collection Containers


Application in three easy steps



1

Prepare the Container


Use a clean, dry urine collection container for sample collection.

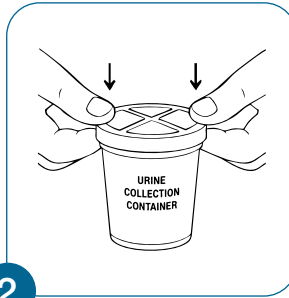

2

Collect the Sample

Collect urine directly into the container without exceeding the 8 oz (240 mL) capacity.

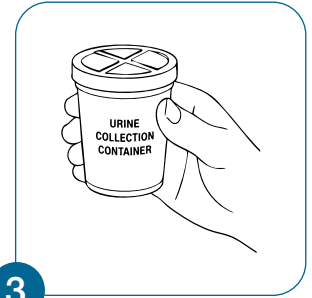

3

Secure the Lid

Place the plastic lid on the container to prevent leakage during transport.

Single-Use Design:

Single use. Applies to Medzah Urine Collection Containers, Item Number: CP10221.

Contents:

280 gsm paper, 18 gsm plastic liner, and 2 gsm plastic lid.

Caution:

For urine specimen collection only. Do not use if the container or lid is damaged or compromised. Do not reuse the container after use. Follow a healthcare professional's directions.

Storage and Handling:

Store in a clean, dry area away from excessive heat, moisture, and direct sunlight. Handle the container carefully to prevent damage. Keep the lid securely closed when transporting a collected specimen.

Important Use Information:

Use a clean container for each specimen collection. Collect the urine specimen directly into the container and securely place the lid after collection. Ensure the container is properly identified and transported according to applicable laboratory procedures.

