


Wrist Brace


Application in three easy steps

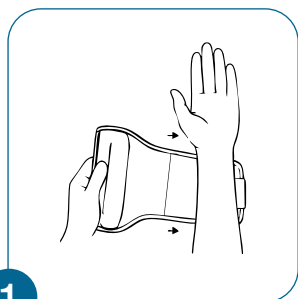

1

Position the Brace


Place the brace around the wrist and adjust the expandable panel for a comfortable fit.

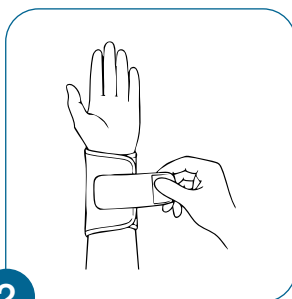

2

Secure the Strap

Fasten the auto-lock strap securely for a firm and comfortable fit.

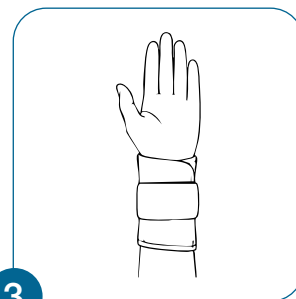

3

Check the Fit

Ensure the brace provides firm wrist support, allowing unrestricted movement.

Reusable Design:

Repeated use. Applies to Medzah Wrist Brace, Item Number: 3000-10, 3000-11, 3000-12, 3000-20, 3000-21, 3000-22.

Contents:

Polyester, nylon, aluminum, polypropylene, steel, and hook-and-loop components.

Caution:

Do not use if the brace or strap is damaged. Ensure the brace is securely fitted without excessive pressure or discomfort. Adjust the brace as needed for a comfortable fit and follow a healthcare professional's directions.

Storage and Handling:

Store at room temperature in a clean, dry location. Keep away from excessive moisture, heat, and sharp objects. Inspect the brace before each use.

Important Use Information:

Select the correct size and side (left or right) before use. Inspect the brace and strap before use. Adjust the expandable panel to conform to the wrist and secure the auto-lock strap for a comfortable, stable fit. Ensure hand and thumb movement remains unrestricted.

