
Disposable 3-Ply Face Masks

Use in three easy steps

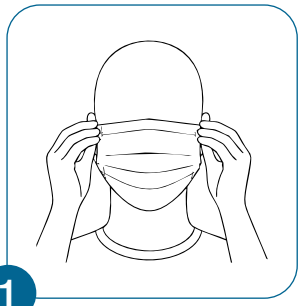

1

Position the Mask

Place the mask over your face, covering your nose, mouth, and chin.



2

Secure the Earloops

Place the earloops around both ears and adjust for a comfortable fit.



3

Check the Fit

Ensure the mask stays in place and fully covers your nose, mouth, and chin.

Single-Use Design:
Single use. Applies to Medzah Disposable 3-Ply Face Masks, Item Number: MDZH-FM01.

Contents:
Disposable 3-ply mask with 25 gsm non-woven and melt-blown layers and elastic earloops. Latex-free. ≥95% bacterial and particulate filtration.

Caution:
Inspect the mask before use and do not use if it is damaged, torn, or contaminated. Ensure the mask is properly positioned and completely covers the nose, mouth, and chin.

Storage and Handling:
Store masks in a clean, dry location away from moisture, excessive heat, and contamination. Keep the masks in their original packaging until ready for use and handle with clean hands.

Important Use Information:
This mask is designed for single use in healthcare and clinical environments. Replace the mask when it becomes damp, soiled, damaged, or difficult to breathe through. Meets ASTM Level 1 requirements.

Scan for product details

medzah.com